



Težina uređaja

1.9kg

Dimenzije uređaja

188x173x163.5mm

Veličina ekrana

127x98mm (5-inch)

Kapacitet ovlaživača

300ml ±20ml

Razina buke

25dB (A)

Daljinski nadzor terapije

Podržava

Automatsko online ažuriranje sustava

OTA

Auto Start/Stop


Podržava

Rampa

0-45min, Auto

Povezivost

Bluetooth/WiFi

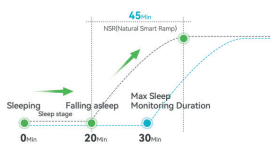

### Pametna prirodna rampa

- Automatsko i precizno prepoznavanje početka sna
- Terapija započinje nižim tlakom koji se postupno i ugodno povećava do zadane vrijednosti
- Maksimalno vrijeme rampe iznosi do 45 minuta

### Dvostruki sustav ovlaživanja

- Neovisni dvoslojni dizajn s učinkovitošću ovlaživanja od 96 %
- Automatska kontrola vlažnosti prilagođava se uvjetima okoline za maksimalnu udobnost korisnika

### Prikaz terapijskih krivulja u stvarnom vremenu

- Trostruki prikaz terapijskih krivulja s podacima o protoku, tlaku i tidalnom volumenu
- Podrška za analizu krivulja uz terapijske podatke u stvarnom vremenu
- Detaljna analiza krivulja uz mogućnost povećanja do 8x

### Inovativan dizajn


- Kompaktne dimenzije omogućuju maksimalnu iskoristivost prostora
- Spremnik za vodu jednostavan za čišćenje s magnetskim zatvaranjem za lakše održavanje
- Upravljanje pritiskom na tipku uz znanstveno optimiziranu konfiguraciju i stabilne performanse

### Automatska kompenzacija curenja zraka

- Detekcija i kompenzacija curenja zraka u stvarnom vremenu
- Održava stabilan tlak u dišnim putovima
- Potpuno automatski rad uređaja

### Sinkronizacija prirodnog disanja

- Brza reakcija uređaja smanjuje otpor pri disanju
- Višerazinske prilagodbe za personaliziranu terapiju
- Povećana udobnost poboljšava suradljivost i pridržavanje terapiji

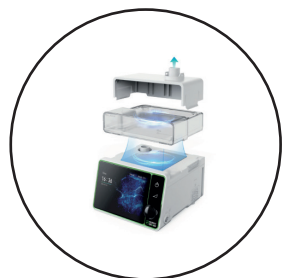

## Potpuno dodirni zaslon (5")

- Zaslon visoke rezolucije s intuitivnim upravljanjem
- Potpuni pregled podataka na prvi pogled
- Upravljanje putem dodira za jednostavnu navigaciju

## Ultra tih rad – 25 dB(A)

- Napredni antivibracijski dizajn s višeslojnim smanjenjem buke
- Inteligentna tehnologija kontrole buke

## OTA daljinska ažuriranja

- Automatska ažuriranja jednim klikom za jednostavnije korištenje
- Pristup najnovijim verzijama softvera u bilo kojem trenutku
- Poboľjšano korisničko iskustvo kroz kontinuirani razvoj softvera.

## Potpuna pametna povezanost

- Integracija više uređaja putem opreme, mobilnih uređaja, računala i clouda
- Platforma SunnyCenter™ i aplikacija SunnyCare™ omogućuju sveobuhvatno upravljanje respiratornom terapijom i praćenjem spavanja putem IoT tehnologije

Model Comparison Table

Product Series		CPAP (Green)	BIPAP (Blue)		NIV (Golden)			
Product Model		NatrSleep 20A	NatrSleep 20S_Auto	NatrSleep 25S_Auto	NatrRes 30ST+	NatrRes 25Vp	NatrRes SVAP Pro	NatrRes 30 Premium
Comfort Setting	Pressure Range	4–20 cmH <sub>2</sub> O	2–20 cmH <sub>2</sub> O	2–25 cmH <sub>2</sub> O	2–30 cmH <sub>2</sub> O	2–25 cmH <sub>2</sub> O	2–30 cmH <sub>2</sub> O	2–30 cmH <sub>2</sub> O
	Inhale Time(s)	/	0.3–4s	0.3–4s	0.3–4s	0.3–4s	0.3–4s	0.3–4s
	Trigger Sensitivity	/	5 levels	5 levels	5 levels	5 levels	5 levels	5 levels
	Humidification	0–6 levels/Auto	0–6 levels/Auto	0–6 levels/Auto	0–6 levels/Auto	0–6 levels/Auto	0–6 levels/Auto	0–6 levels/Auto
Operation Mode	CPAP	●	●	●	●	●	●	●
	APAP	●						
	S		●	●	●		●	●
	S-Auto		●	●	●			
	T						●	●
	S/T					●		●
	ASV						●	●
	ASV Pro						●	●
	SPC							●
SVAPS						●	●	
Data Management	WiFi	●	●	●	●	●	●	●
	Bluetooth	●	●	●	●	●	●	●
	SD Card	●	●	●	●	●	●	●
	Mobile APP	●	●	●	●	●	●	●
	Cloud Platform	●	●	●	●	●	●	●
Treatment Curve		AHI Trend		P95 Trend		Intelligent Alarm		
Support		Support		Support		Support		

# NatrSync PSGmate

Potpuno automatski titracijski sustav

PRECIZNOST SPOJENA S VRHUNSKIM PERFORMANSAMA

Podržavajući 10 načina disanja i 13 parametara praćenja, uređaj donosi modularni dizajn za brzo i jednostavno postavljanje. Spaja jednostavnost i naprednu funkcionalnost kroz integrirani dizajn i univerzalnu kompatibilnost sa svim PSG sustavima.

## Kompaktan i snažan

- Kompaktan „plug-and-use“ dizajn spreman za jednostavno korištenje
- Modularni dizajn uklanja potrebu za dodatnim naponskim vodovima
- 16 integriranih priključaka kompatibilnih sa svim PSG uređajima
- Jednolinijsko povezivanje za jednostavnije kabliranje
- Wi-Fi povezanost za neometanu mobilnost unutar bolničkih odjela

## Pametni sustav upravljanja titracijom

- Instalacija softvera jednim klikom za automatsku konfiguraciju
- Fleksibilni načini rada za titraciju jednog ili više pacijenata i dodatne mogućnosti prilagodbe
- Prilagodljivo upravljanje podacima pacijenata uz automatiziranu analizu parametara praćenja
- Podržava više terminalnih uređaja: računala, tablete, pametne telefone i daljinsko upravljanje

